




Test Catalog

Diagnostic. Prognostic. Predictive. Predisposition.



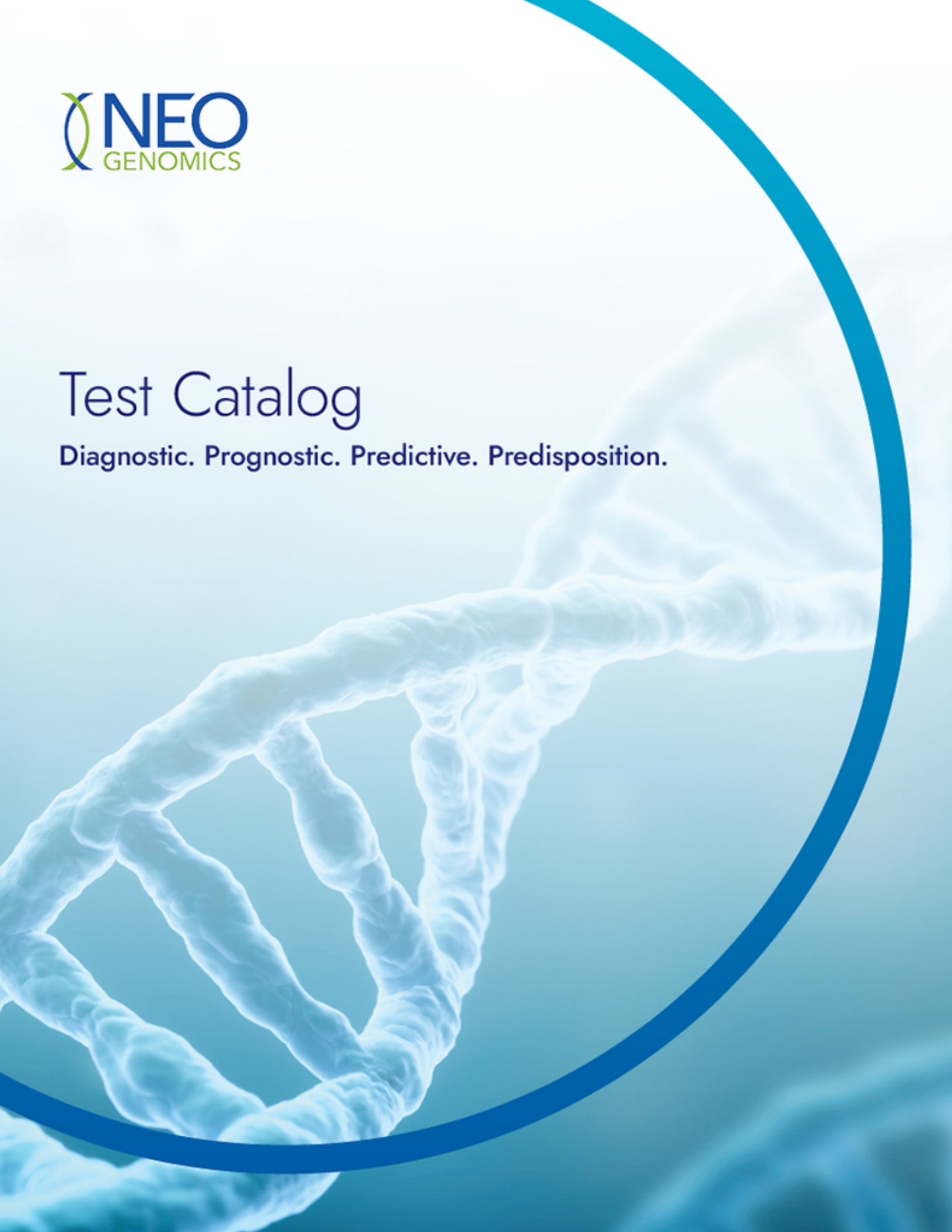

NeoTYPE® Head & Neck Tumor Profile

Alternative Name

Head & Neck Tumor Profile

Methodology

Molecular

Test Description

The NeoTYPE Head & Neck Tumor Profile analyzes 34 biomarkers through a combination of next-generation sequencing (NGS), other molecular methods, FISH, and IHC as listed below. Test orders include summary interpretation of all results that help guide treatment decisions. If Pan-TRK IHC is expressed or equivocal, reflex to either NTRK NGS Fusion Panel (Default) or NTRK 1-3 FISH Panel will be added. A microsatellite instability (MSI) NGS result of "indeterminate" will create a reflex to MSI by PCR as long as the tumor percentage is ≥40% and paired normal tissue is available.

- NGS (27 genes + 2 biomarkers): AKT1, ATM, BRAF, CDKN2A, CTNNB1, EGFR, ERBB2, ERBB4, FBXW7, FGFR1, FGFR2, FGFR3, HRAS, IDH1, IDH2, KRAS, MET, Microsatellite Instability (MSI), NFE2L2, NOTCH1, NRAS, PIK3CA, PTEN, RB1, SMO, SRC, TERT Promoter, TP53, Tumor Mutation Burden (TMB)
- Other Methodology (1 biomarker): HPV RNA ISH 16/18 (High Risk)
- FISH (3 FISH): MET, PTEN, RET (tech-only available)
- IHC (2 biomarkers): PD-L1 22C3, Pan-TRK (tech-only available for PD-L1)

Clinical Significance

The NeoTYPE Head and Neck Tumor Profile characterizes head and neck tumors of any histological subtype for the most significant genetic changes relevant to therapy decisions, prognosis, and clinical research. It is appropriate for patients with newly-diagnosed or recurrent disease, and for patients with resistant disease to explore options in clinical trials.

Specimen Requirements

- **FFPE tissue:** Paraffin block preferred. Please use 10% buffered formalin fixative. Do not use zinc fixatives.

Storage & Transportation

Use cold pack for transport, making sure cold pack is not in direct contact with specimen.

CPT Code(s)*

81445x1, 88377x3, 88360x1, 88342x1, 88365x1; add 81194x1 if reflexed to NTRK NGS Fusion Panel (default) or 88374x3 automated (88377x3 manual) if reflexed to NTRK 1-3 FISH Panel

Medicare MoIDX CPT Code(s)*

81479x1, 88377x3, 88360x1, 88342x1, 88365x1; add 81194x1 if reflexed to NTRK NGS Fusion Panel (default) or 88374x3 automated (88377x3 manual) if reflexed to NTRK 1-3 FISH Panel

New York Approved

Yes

Level of Service

Global

Turnaround Time

14 days

*The CPT codes provided with our test descriptions are based on AMA guidelines and are for informational purposes only. Correct CPT coding is the sole responsibility of the billing party.

Please direct any questions regarding coding to the payor being billed.

NeoGenomics Laboratories is a specialized oncology reference laboratory providing the latest technologies, testing partnership opportunities, and interactive education to the oncology and pathology communities. We offer the complete spectrum of diagnostic services in molecular testing, FISH, cytogenetics, flow cytometry, and immunohistochemistry through our nation-wide network of CAP-accredited, CLIA-certified laboratories.

Committed to research as the means to improve patient care, we provide Pharma Services for pharmaceutical companies, in vitro diagnostic manufacturers, and academic scientist-clinicians. We promote joint publications with our client physicians. NeoGenomics welcomes your inquiries for collaborations. Please contact us for more information.

*The CPT codes provided with our test descriptions are based on AMA guidelines and are for informational purposes only. Correct CPT coding is the sole responsibility of the billing party.

Please direct any questions regarding coding to the payor being billed.



9490 NeoGenomics Way
Fort Myers, FL 33912
Phone: 239.768.0600/ Fax: 239.690.4237
neogenomics.com
© 2022 NeoGenomics Laboratories, Inc. All Rights Reserved.
All other trademarks are the property of their respective owners
Rev. 120922