




# Test Catalog

Diagnostic. Prognostic. Predictive. Predisposition.



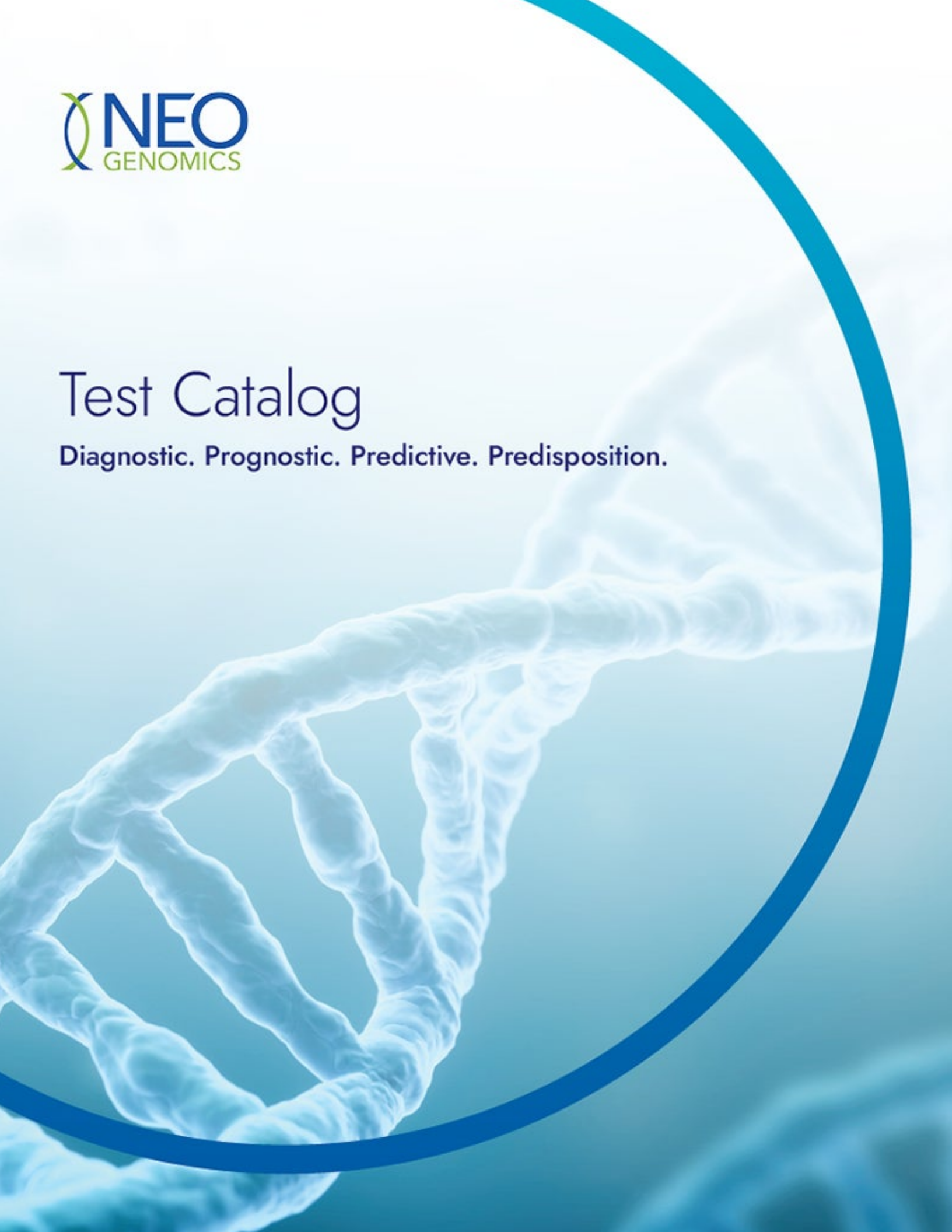

## ALL Adult FISH Panel

### Alternative Name

Acute lymphoblastic leukemia

### Methodology

FISH

### Test Description

**Probes:** TCF3/PBX1 (E2A/PBX1) t(1;19) | Trisomy or Tetrasomy 4, 6, 10, 17 (Cen 4, Cen 6, Cen 10, Cen 17) | MYC (8q24) | BCR/ABL1/ASS1 t(9;22) | MLL (11q23) | IgH (14q32) |

**Disease(s):** Acute lymphoblastic (lymphocytic) leukemia (B-cell ALL), B lymphoblastic lymphoma (LBL), adult  
Probes may be ordered separately except Centromeres 4 and 17 are paired, and Centromeres 6 and 10 are paired.

**Note:** STAT processing is available by request for BCR-ABL1. Note STAT along with MD contact name and phone number to receive STAT results.

**Note:** CDKN2A (p16) Deletion FISH is also available and may be ordered separately. See details [here](#).

### Clinical Significance

The ALL Adult FISH Panel is used for the detection of recurrent chromosome abnormalities observed in adults with ALL of B-cell lineage and B lymphoblastic lymphoma (LBL). Identification of specific abnormalities helps predict disease aggressiveness and response to therapy. This panel differs from the ALL Pediatric FISH Panel in that this panel excludes probes for ETV6/RUNX1 t(12;21).

### Specimen Requirements

- **Bone Marrow Aspirate:** 1-2 mL sodium heparin tube. EDTA tube is acceptable.
- **Peripheral Blood:** 2-5 mL sodium heparin tube. EDTA tube is acceptable.
- **Fresh, Unfixed Tissue:** Tissue in RPMI.
- **Fluids:** Equal parts RPMI to specimen volume.
- **Paraffin Block or Cut Slides:** Not available.
- **Note:** Please exclude biopsy needles, blades, and other foreign objects from transport tubes. These can compromise specimen viability and yield, and create hazards for employees.

### Storage & Transportation

Use cold pack for transport, making sure cold pack is not in direct contact with specimen.

### CPT Code(s)\*

88374x7 automated. Codes may differ if manual analysis is performed.

### New York Approved

Yes

### Level of Service

Global, Technical

## Turnaround Time

3-5 days. STAT results for BCR-ABL1, when requested, are reported 12-24 hours from receipt in the NeoGenomics laboratory.

\*The CPT codes provided with our test descriptions are based on AMA guidelines and are for informational purposes only. Correct CPT coding is the sole responsibility of the billing party.

Please direct any questions regarding coding to the payor being billed.

NeoGenomics Laboratories is a specialized oncology reference laboratory providing the latest technologies, testing partnership opportunities, and interactive education to the oncology and pathology communities. We offer the complete spectrum of diagnostic services in molecular testing, FISH, cytogenetics, flow cytometry, and immunohistochemistry through our nation-wide network of CAP-accredited, CLIA-certified laboratories.

Committed to research as the means to improve patient care, we provide Pharma Services for pharmaceutical companies, in vitro diagnostic manufacturers, and academic scientist-clinicians. We promote joint publications with our client physicians. NeoGenomics welcomes your inquiries for collaborations. Please contact us for more information.

\*The CPT codes provided with our test descriptions are based on AMA guidelines and are for informational purposes only. Correct CPT coding is the sole responsibility of the billing party.

Please direct any questions regarding coding to the payor being billed.



9490 NeoGenomics Way  
Fort Myers, FL 33912  
Phone: 239.768.0600/ Fax: 239.690.4237  
neogenomics.com  
© 2023 NeoGenomics Laboratories, Inc. All Rights Reserved.  
All other trademarks are the property of their respective owners  
Rev. 012823